


Y27632

Product Description

Y27632 is a cell-permeable small molecule Rho-associated kinase (ROCK) inhibitor ^[1]. Y27632 has been found to prevent apoptosis as well as enhance the survival and cloning efficiency of dissociated human embryonic stem (ES) cells without affecting their self-renewal properties or pluripotency ^[2]. This molecule has also been shown to enhance survival during the transplantation of ES cell-derived neural precursors ^[3]. Y27632 in combination with Pifithrin- μ significantly improves cell recovery after cryopreservation ^[4].

Catalog No.:	MD-0025
Quantity:	2mg (10mM in 624.4 μ l DMSO)
Alternate Name:	(1R,4r)-4-((R)-1-aminoethyl)-N-(pyridine-4-yl)cyclohexanecarboxamide dihydrochloride
CAS Number:	146986-50-7
Chemical Formula:	C ₁₄ H ₂₃ Cl ₂ N ₃ O
Molecular Weight:	320.3 g/mol
Purity:	>Greater than 98% by HPLC analysis

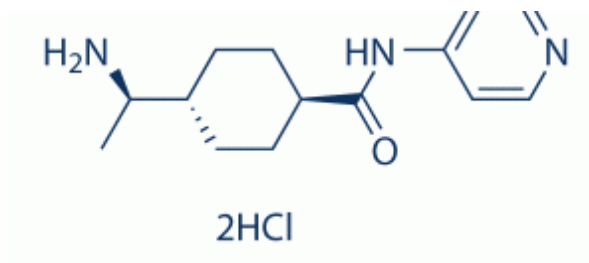

Structure:

Applications: MAINTENANCE AND SELF-RENEWAL

- Enhances survival of human embryonic stem (ES) cells when they are dissociated to single cells by preventing dissociation-induced apoptosis (anoikis), thus increasing their cloning efficiency.
- Improves embryoid body formation using forced-aggregation protocols.
- Increases the survival of cryopreserved single human ES cells after thawing.
- Blocks apoptosis of mouse ES-derived neural precursors after dissociation and transplantation.

REPROGRAMMING

- Direct lineage reprogramming of fibroblasts to mature neurons, in combination with CHIR99021, RepSox, Forskolin, SP600125, Gö6983 and Valproic Acid.

DIFFERENTIATION

- Improves survival of human ES cell monolayers at the initiation of differentiation protocols

Product Use: Y27632 is for research use only. It is not approved for human or animal use, or for application in *in vitro* diagnostic procedures.

Storage & stability: Store aliquots at -20°C. Stock solutions are stable for 6 months when stored as directed.

References

[1] Ishizaki, T., Uehata, M., Tamechika, I., Keel, J., Nonomura, K., Maekawa, M., and Narumiya, S. (2000) Pharmacological properties of Y-27632, a specific inhibitor of rho-associated kinases. *Mol Pharmacol* 57: 976-983.

[2] Watanabe, K., Ueno, M., Kamiya, D., Nishiyama, A., Matsumura, M., Wataya, T., Takahashi, J.B., Nishikawa, S., Nishikawa, S., Muguruma, K., and Sasai, Y. (2007) A ROCK inhibitor permits survival of dissociated human embryonic stem cells. *Nat Biotechnol* 25: 681-686.

[3] Koyanagi, M., Takahashi, J., Arakawa, Y., Doi, D., Fukuda, H., Hayashi, H., Narumiya, S., and Hashimoto, N. (2008) Inhibition of the Rho/ROCK pathway reduces apoptosis during transplantation of embryonic stem cell-derived neural precursors. *J. Neurosci Res.* 86: 270-280.

[4] Xu, X., Cowley, S., Flaim, C.J., James, W., Seymour, L.W., Cui, Z. (2010) Enhancement of cell recovery for dissociated human embryonic stem cells after cryopreservation. *Biotechnol Prog.* 26: 781-788.

Disclaimers

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures.

© iXCells Biotechnologies 2015. All rights reserved.